

HRT-7000 ADVANCED REFRACTION TABLE





HRT-7000

ADVANCED REFRACTION TABLE

Another frontier step in refraction table system with unique handle mechanism: more comfortable and convenient HRT-7000



HANDLE DESIGN

Ergonomic design with unique handle mechanism to facilitate the movement of upper table with extreme ease



POWER ARM

45 degree soft movement and up & down multi sequence adjustment will enable to accommodate any patient comfortably.




UP & DOWN FUNCTION

Up & Down feature is added in order to provide convenience for measuring position of patient's eyes.



KEYPAD

Visible location of keypads enables you to control the power arm and chair precisely with ease.

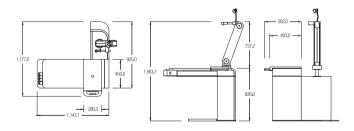

Designs and details are subject to change without prior notice for their improvements.



Huvitz Building 689-3, Geumjeong-dong Gunpo-si, Gyeonggi-do, 435-862, Korea

Tel: +82-31-442-8868 Fax: +82-31-477-8617 http://www.huvitz.com

SPECIFICATIONS

Dimensions	1,147(W) x 1,177(D) x 1,564(H)mm
Table Size	1,147(W) x 450(D) x 30(T)mm
Table Rotating Angle	90°
Head Stroke	400mm(R-L)
Power Arm Up-Down Range	±100mm
Chair Up-Down Range	150mm
Lamp	LED Lamp 3W
Power Supply	AC 110/220V, 50/60Hz
Power Consumption	420W (Max)

Distributed by