
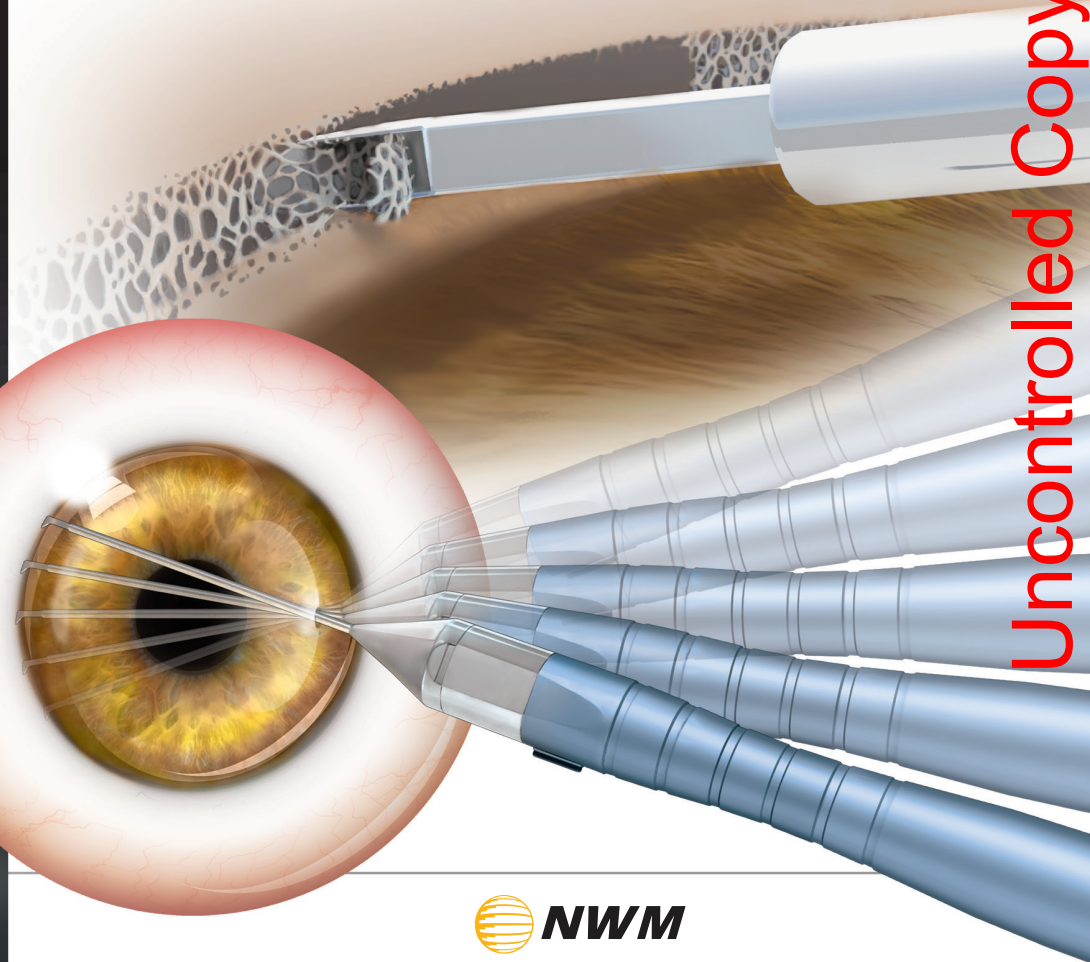


KAHOOK

DUAL BLADE

Features:

- Single use, ophthalmic blade
- Utilizes ab interno approach through 1.2mm clear corneal incision
- Unique blades designed to fit in canal of Schlemm
- Precisely placed dual blades to excise trabecular meshwork and minimize residual leaflets
- Maintains natural physiologic outflow pathways

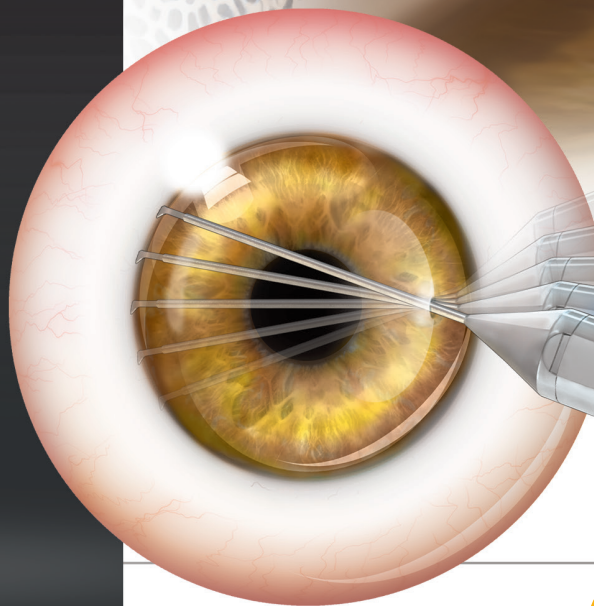

909.466.4304

10763 Edison Court
Rancho Cucamonga, CA 91730, USA


Uncontrolled Copy If Printed



DUAL BLADE



1. Pointed tip easily pierces the trabecular meshwork
2. Ramp elevates and stretches trabecular meshwork
3. Dual blades excise strip of trabecular meshwork
4. Foot plate prevents damage to anterior wall of canal and facilitates smooth motion

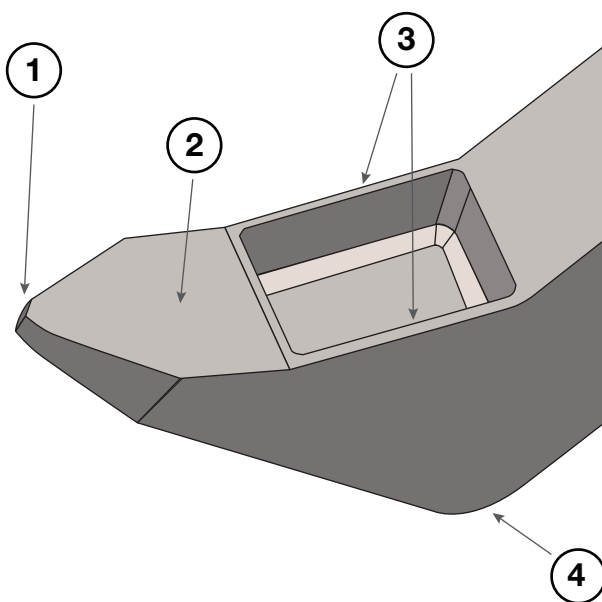

**Get Certified
for KDB surgery @**
KDBcert.com



New World Medical, inc.
909.466.4304

www.NewWorldMedical.com

10763 Edison Court
Rancho Cucamonga, CA 91730

Part # 50-0076 Rev. New 10/15

Uncontrolled Copy If Printed