





CANABRAVA RING

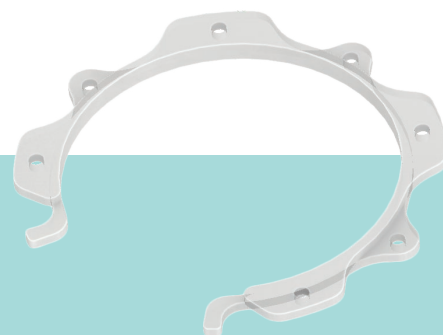

Iris Expander Ring

Brand name	CANABRAVA RING
Reference	CANA - RING
Description and Indications	<p>An iris expander ring designed to safely expand the iris in order to simplify phacoemulsification, as well as surgeries in which the pupil does not dilate properly using pharmaceuticals. It provides a pupillary opening of 6.3mm and keeps dilation stable during surgery, improving the surgeon's access and field of vision.</p> <ul style="list-style-type: none">• 300° semicircular ring with a 60° opening made from a single piece.• Does not require an injector.• Easy to adapt and remove.• The ring can be inserted in incisions smaller than 2 mm.• Easy to insert in pupils smaller than 4.5 mm and in shallow chambers.• Optimum visibility.• Disposable device that cannot be re-sterilized. <ul style="list-style-type: none">• Myosis.• Coloboma.• Traumatic sector iridectomy.• Intraoperative floppy iris syndrome.• Iris prolapse.• Small pupils.• Rigid iris tissue due to a pseudoexfoliation.• Vitrectomy with small pupils.
Material	Medical grade Polymethyl Methacrylate (PMMA).
Sterilization method	Ethylene Oxide (ETO).
Packaging	Sterile, single unit.
Shelf-life	4.5 years

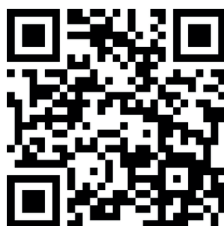

REV: 05. 07/2022