





CANABRAVA RING

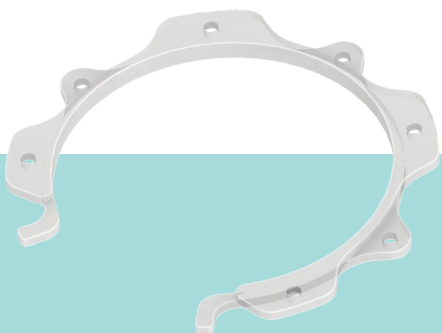

Anillo Extensor de Iris


Nombre de marca	CANABRAVA RING
Referencia	CANA - RING
Descripción e indicaciones	El anillo extensor de iris está diseñado para proporcionar una dilatación pupilar estable y sostenida mediante la extensión mecánica del iris, con el fin de facilitar el campo de visión del cirujano y mejorar el acceso a la cámara posterior durante la cirugía oftalmológica. • Anillo semicircular de 300º con apertura de 60º fabricado de una sola pieza. • No precisa inyector. • Sencilla adaptación y extracción. • Fácil colocación en pupilas menores de 4,5 mm y en cámaras estrechas. • Óptima visualización. • Producto de un solo uso Indicaciones: • Miosis o pupilas pequeñas. • Coloboma de iris. • iridectomía sectorial trumática. • Síndrome de iris flácido intraoperativo. • Prolapso de iris. • Tejido del iris rígido debido a una pseudoexfoliación. • Vitrectomía con pupilas pequeñas.
Material	Polimetilmetacrilato de grado médico (PMMA).
Método de esterilización	Óxido de etileno (ETO).
Presentación	Envasado unitario estéril.
Caducidad	4,5 años

CE 0318

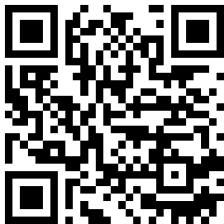



CANABRAVA RING

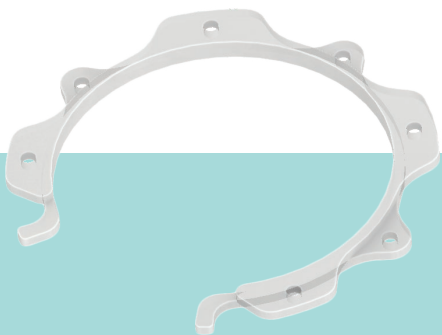

Iris Expander Ring

Brand name	CANABRAVA RING
Reference	CANA - RING
Description and Indications	The Iris Extender Ring is designed to provide stable and sustained pupillary dilation by mechanically extending the iris to facilitate the surgeon's field of vision and improve access to the posterior chamber during ophthalmic surgery. • 300° semicircular ring with a 60° opening made from a single piece. • Does not require an injector. • Easy to adapt and remove. • Easy to insert in pupils smaller than 4.5 mm and in shallow chambers. • Optimum visibility. • Single-use product Indications: • Myosis or small pupils. • Iris coloboma. • Traumatic sectorial iridectomy. • Intraoperative floppy iris syndrome. • Prolapse of iris. • Tissue of the rigid iris due to a pseudoexfoliation. • Vitrectomy with small pupils.
Material	Medical grade Polymethyl Methacrylate (PMMA).
Sterilization method	Ethylene Oxide (ETO).
Packaging	Sterile, single unit.
Shelf-life	4.5 years

CE 0318

